



Thermo Scientific NanoDrop Products

Instruments for rapid microvolume analysis of DNA, RNA and protein samples

NanoDrop—trusted by scientists worldwide

Rely on fast, accurate quantification of DNA, RNA and protein samples using only 1–2 μL with Thermo Scientific™ NanoDrop™ microvolume instruments. No dilutions needed even for highly concentrated samples with pioneering sample-retention technology* using optical measurement pedestals. For over 15 years and with over 40,000 citations, NanoDrop instruments have been helping scientists around the world do their best work. With preconfigured methods for common life science applications NanoDrop instruments make a novice perform like an expert. It's as simple as pipette, measure, know!



NanoDrop One/One^C UV-Vis Spectrophotometer

Intelligent analysis, streamlined workflows

- Contaminant identification and corrected concentrations with Thermo Scientific™ Acclaro™ Sample Intelligence technology
- Touchscreen interface with guided troubleshooting and pre-programmed methods for life sciences
- Modern connectivity allows seamless data transfer to PC or network via Wi-Fi, USB, or Ethernet. Export your data to Thermo Fisher Connect and access it anytime, anywhere from any device!

NanoDrop products comparison—
**choose the instrument
that's right for you**

Instrument	Measures 1–2 μL sample	Pre-programmed methods for life sciences ¹	Full-spectral data	NUCLEIC ACIDS	
				A260	A260/A280
NanoDrop One/One ^C	UV-Vis	✓	✓	✓	✓
NanoDrop 8000	UV-Vis	✓	✓	✓	✓
NanoDrop Lite	UV-Vis	✓		✓	✓
NanoDrop 3300	Fluorescence	✓	✓		

1. Bradford, BCA, Lowry, Pierce 660 2. NanoDrop One and NanoDrop 8000 use absorbance to measure fluorescently-labelled nucleic acids and proteins



NanoDrop 8000 UV-Vis Spectrophotometer

Higher throughput, full-spectrum analysis


- Measure 1 to 8 samples simultaneously on a linear array of pedestals
- Simply pipette directly from tubes or 96-well plates. Analyze 96 samples in <6 minutes
- Preconfigured methods for common applications such as Nucleic Acids, Protein A280, Microarray, Proteins and Labels, Pierce 660, Bradford, BCA, Lowry and more



NanoDrop Lite UV-Vis Spectrophotometer

Simple analysis, compact delivery

- Dedicated 260 nm and 280 nm wavelengths to measure nucleic acid and protein concentrations
- Palm-sized, portable instrument requires no PC control
- Optional docking printer prints cryogenic labels

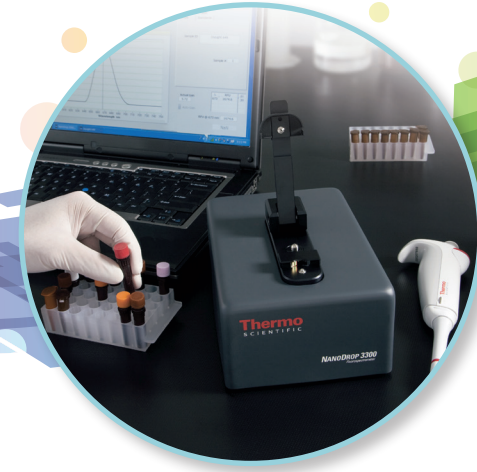

NanoDrop 3300 Fluorospectrometer

Microvolume fluorescence analysis


- Wide excitation range without filter changes or monochromators
- Diverse range of measurements from routine (e.g., Thermo Scientific™ RiboGreen™ fluor, Thermo Scientific™ PicoGreen™ fluor, GFP) to advanced fluorescence analysis (e.g., nanoparticles, quantum dots, FRET)

NUCLEIC ACIDS SAMPLES			PROTEIN AND PEPTIDE SAMPLES					Custom methods editor	High-res touchscreen interface	Requires computer to operate	Can evaluate FLR-labeled samples ²	Advanced connectivity ³	Auto-Measure capability
A280	A260/A230	Acclaro contaminant ID	A280	A260/A280	A205	Colorimetric assays ¹	Acclaro contaminant ID						
	✓	✓	✓	✓	✓	✓	✓	✓		✓	✓	✓	
	✓		✓	✓		✓		✓	✓		✓		
	✓		✓	✓						✓			
								✓		✓			

¹ d proteins ³ USB, Wi-Fi and Ethernet connections to PC or network, print from local control

²US Patents 6628382 and 6809826

NanoDrop instrument specifications

Description	NanoDrop One/One [®] Spectrophotometer	NanoDrop Lite Spectrophotometer	NanoDrop 8000 Spectrophotometer	NanoDrop 3300 Fluorospectrometer
Instrument Control	Built-in touchscreen	Local control	Computer software	Computer software
Minimum Sample Size	1 µL	1 µL	1 µL	1 µL
Sample Number	1	1	up to 8	1
Pathlength(s)	0.030 to 1.0 mm auto-ranging	0.5 mm	1.0 to 0.2 mm, auto-ranging	N/A
Light Source(s)	Xenon flash lamp	Light emitting diodes	Xenon flash lamp	Light emitting diodes
Excitation Maxima of LEDs	N/A	N/A	N/A	UV: 365 nm, Blue: 470 nm, White: 460–650 nm
Detector Type	2048-element CMOS linear image sensor	Silicon photodiode	2048-element linear silicon CCD array	2048-element linear silicon CCD array
Wavelength Range	190–850 nm	260 and 280 nm	220–750 nm	400–750 nm
Wavelength Accuracy	±1 nm	N/A	±1 nm	±1 nm
Spectral Resolution	≤1.8 nm (FWHM at Hg 254 nm)	<8.0 nm	3 nm (FWHM at Hg 546.1 nm)	8.0 nm (FWHM at Hg 546.1 nm)
Typical Measurement Repeatability	0.002 A (1.0 mm path) or 1%CV, whichever is greater	0.002*	0.003*	<5% CV (10 nM fluorescein)
Absorbance Accuracy**	3% (at 0.97 A, 302 nm)	3% (at 1.05 Abs at 260 nm)	3% (at 0.74 Abs at 350 nm)	N/A
Absorbance Range (10 mm equivalent)	Pedestal: 0–550A, Cuvette: 0–1.5A	0–30 Abs	0–75 Abs	N/A
Lower Limit of Detection	Pedestal: 2 ng/µL (dsDNA) 0.06 mg/mL (BSA) Cuvette: 0.2 ng/µL (dsDNA) 0.006 mg/mL (BSA)	4 ng/µL (dsDNA) 0.12 mg/mL (BSA)	2.5 ng/µL (dsDNA) 0.15 mg/mL (BSA)	<1 fmol fluorescein
Maximum Concentration	Pedestal: 27,500 ng/µL (dsDNA) 820 mg/mL (BSA)	1,500 ng/µL (dsDNA) 45 mg/mL (BSA)	3,700 ng/µL (dsDNA) 100 mg/mL (BSA)	N/A
Measurement and Data Processing Time	8 seconds	<5 seconds	<20 seconds	2–10 seconds
Footprint	20 × 25.4 cm	16 × 11.5 cm	24 × 32 cm	14 × 20 cm
Weight	3.6 kg	0.8 kg	3.4 kg	1.5 kg
Sample Pedestal	303 stainless steel and quartz fiber	303 stainless steel and quartz fiber	303 stainless steel and quartz fiber	303 stainless steel and quartz fiber
Cuvette Position	Optional (with stirring)	N/A	N/A	N/A
Operating Voltage	12 V (DC)	6 V (DC)	12 V (DC)	5 V (DC)
Operating Power Consumption	12–18 W	18 W	30 W	2 W
Standby Power Consumption	5 W	<2.5 W	6 W	1 W
Software Compatibility	Standalone control. PC Viewer: Windows [®] 7 and 10 Professional (64 bit)	Standalone control with data export via USB Flash Drive	Windows [®] 7 and 10 Professional (64 bit)	Windows [®] 7 and 10 Professional (64 bit)
Advanced Connectivity	NanoDrop One Cloud Application*** Powered by  ThermoFisher Cloud	N/A	N/A	N/A

* SD of 10 individual measurements at 0.74 Abs

** Absorbance expressed as Abs/mm measured at 25 °C

*** Visit www.thermofisher.com/connect for details

Find out more at www.thermofisher.com/NanoDrop

ThermoFisher
SCIENTIFIC

© 2017 Thermo Fisher Scientific Inc. All rights reserved. Bluetooth is a registered trademark of Bluetooth SIG, Inc. Windows and Vista are registered trademarks of Microsoft Corporation. All Other trademarks are the property of Thermo Fisher Scientific and its subsidiaries. All NanoDrop instruments are approved to CE and UL/CSA. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details. **BR51530_E 08/17M**